



Certificate of Analysis

GHK Copper 50mg


copper;(2S)-6-amino-2-[[[(2S)-2-[(2-aminoacetyl)amino]-3-(1H-imidazol-5-yl)propanoyl]amino]hexanoate

Compound : GHK-Cu **Client** : Biohacking360 Labs
Lot number : s1-GHK50-742
Analysis date : 2026-03-13
Purity % : 99.67%
Method : HPLC-UV-MS

PubChem CID: 71587328

<https://pubchem.ncbi.nlm.nih.gov/compound/71587328>

High Performance Liquid Chromatography (HPLC) UV – Purity Test

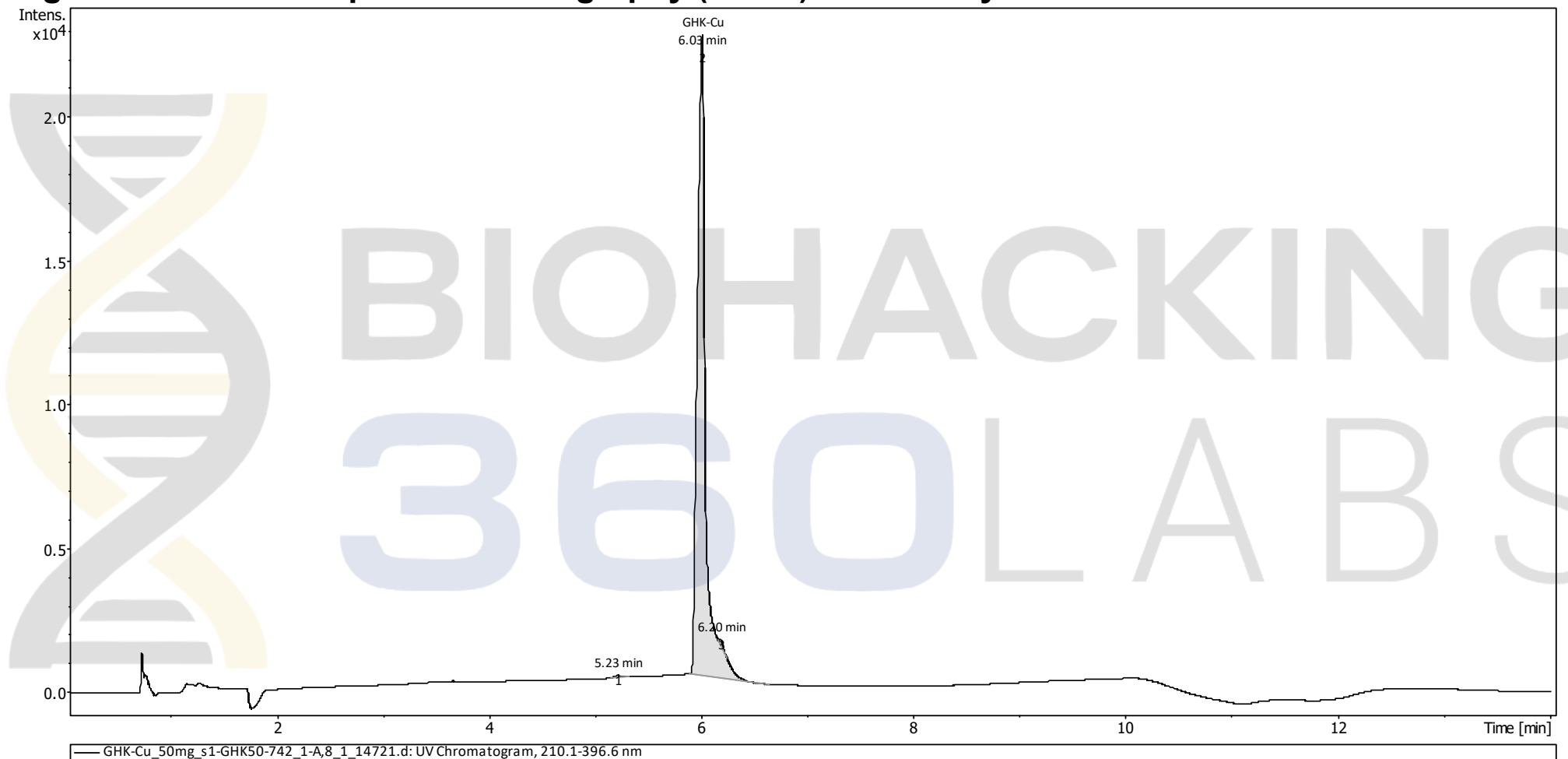

PEAK LIST		Number of detected peaks: 3		
	Time (min)	Area	%Area	
1	5.23	3.48E+02	0.34	
2	6.03	1.02E+05	99.09	GHK Cu
3	6.20	5.91E+02	0.58	GHK

Overall purity : 99.67

The GHK-Copper complex dissociates slightly during HPLC resulting in the small GHK peak immediately after GHK-Cu. This is normal.

Analysis Performed by Ken Pendarvis, ChE Analytical Chemist
 MZ Biolabs
info@biohacking360labs.com

Note: Injectable peptides may contain salts and sugars to aid in solubility and act as pH buffers. These are not normally detected using UV and are not considered impurities.

GHK Copper 50mg

copper;(2S)-6-amino-2-[[[(2S)-2-[(2-aminoacetyl)amino]-3-(1H-imidazol-5-yl)propanoyl]amino]hexanoate

PubChem CID: 73587

<https://pubchem.ncbi.nlm.nih.gov/compound/73587>

Mass Spectrometry (MS) – Identity Test

Identity confirmed using HPLC-MS

Molecular weight calculated using monoisotopic m/z values from mass spectrum

Expected monoisotopic mass : 402.10 Da


Measured monoisotopic mass : 402.12 Da

Molecular weight confirmed

Note : Monoisotopic m/z values are not easily seen in full spectrum view for larger molecules and peptides.

The dominant isotopic peak (base peak) shown in the spectrum below can be used to approximate the average molecular weight frequently reported by vendors and databases as a secondary means of confirmation.

Recorded MS spectrum

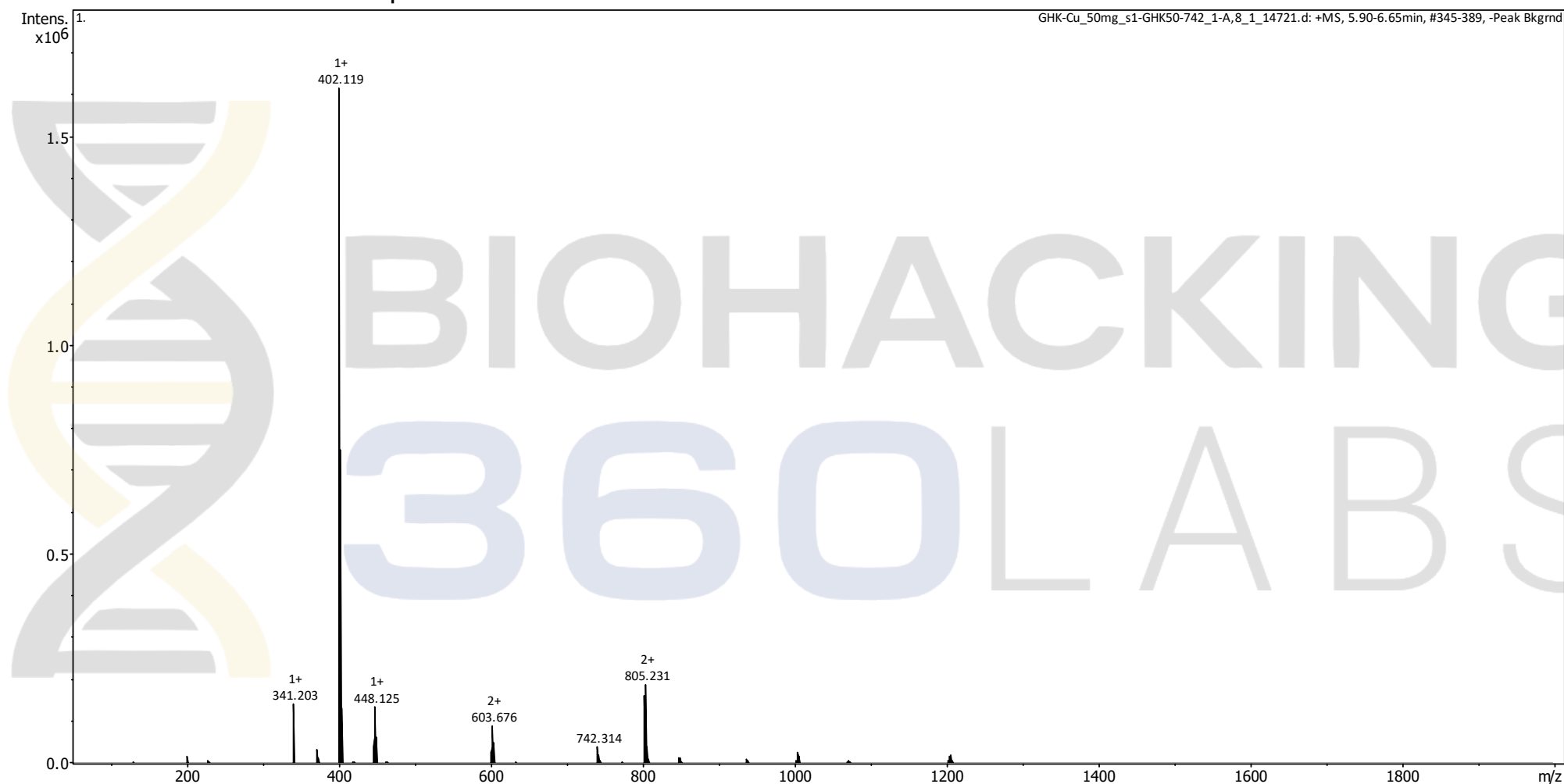

Analysis Performed by Ken
Pendarvis, ChE Analytical
Chemist
MZ Biolabs
info@biohacking360labs.com

